


## SED SPENCER EXTRICATION DEVICE

Extrication/spine immobilization device

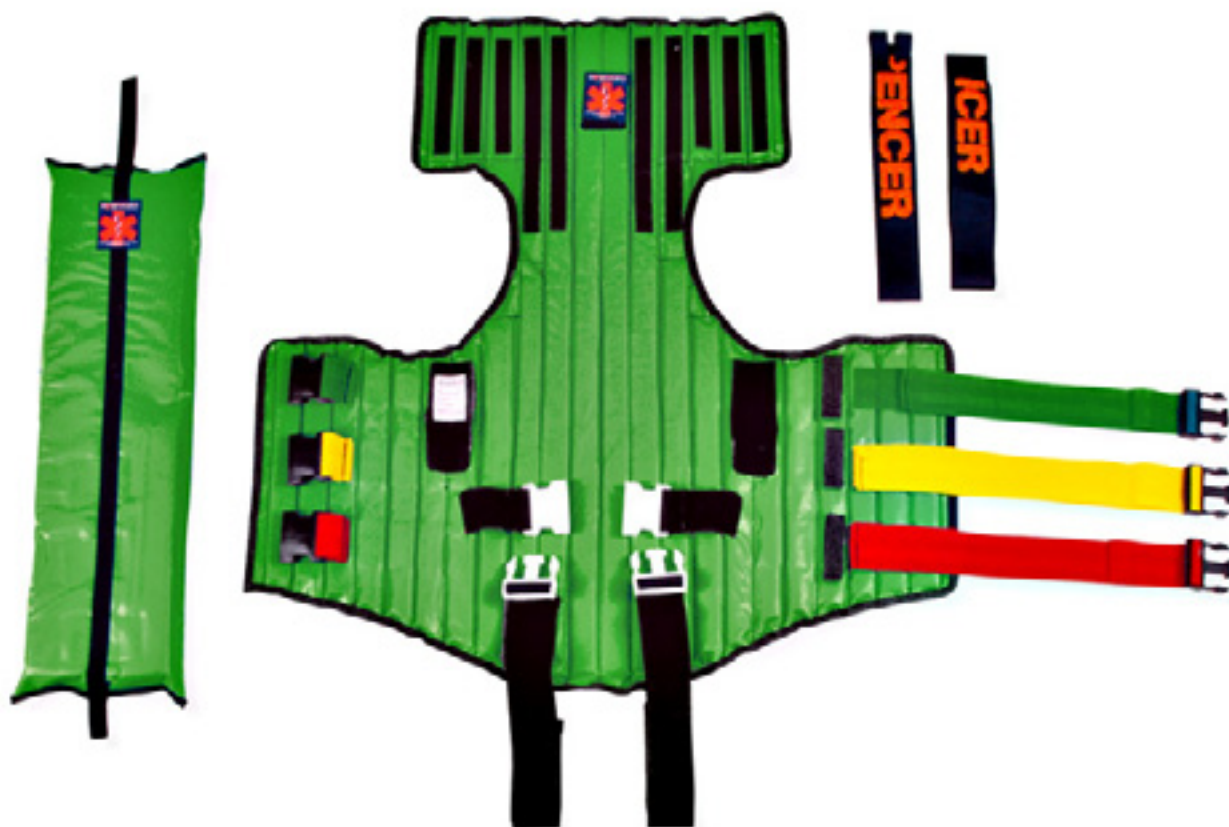

Spencer Extrication Device is an instrument for extrication and spine immobilization, ideal for all emergency situations. Particularly suitable for victims blocked in vehicles, collapsed buildings and disaster scenes. It can also be used with pregnant women, children and babies. It is ideal for the immobilization of the hip and pelvis fractures.

The covering made of robust nylon covered in vinyl limits abrasions and assures total integrity.

The vertical rigidity and the horizontal flexibility are the main characteristics of this device, which allows a rapid immobilization of the spinal column and a secure extrication. The belts and the buckles are coloured for easy recognition. SED is completely X-ray compatible.

It is supplied with

- comfortable and robust transport bag
- nape and chin strap
- cushion to be used as padding for the back of the neck

### TECHNICAL DATA

Dimensions:	860 x 820 mm
Weight with bag and accessories:	3 kg
X-ray compatible:	Yes
Latex free	
NATO stock n°:	6515-15-149-8710

SED Spencer Extrication  
Device SR00111B

